



Hospital Beds and Safety Information for the Patient

A Guide for Patients in the Home

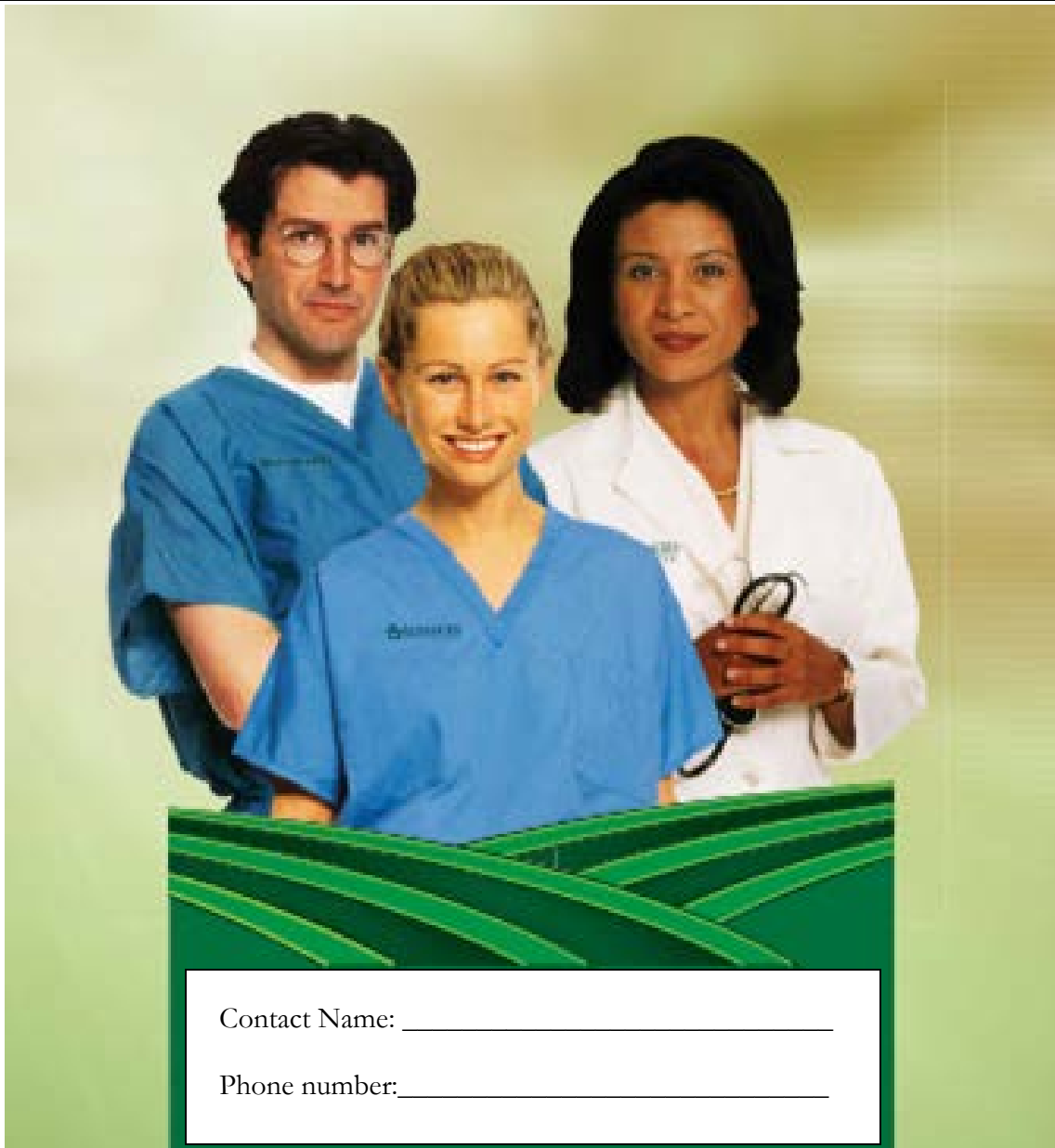




Table of Contents

Chapter 1.....	1
Hospital Bed: Introduction.....	1
Chapter 2.....	2
Instructions for your Hospital Bed	2
What you should do?.....	2
Shaft bed.....	2
Cord bed.....	3
Bed Rails.....	4
Chapter 3.....	5
Special Instructions	5

ICON KEY

 What to do Questions & Answers

Hospital Bed: Introduction

Your doctor has ordered for you to receive a hospital bed. We know you have a choice and are glad you chose Advanced Home Care for your Hospital Bed.





In order to serve you better, we would like to make it easier for you to operate your new hospital bed. Below you will see pictures of several different models of hospital beds Advanced Home Care has to offer. Take a few minutes to get familiar with your hospital bed. We hope this helps you if the power goes out and you need to operate the bed

Instructions for your Hospital Bed




What you should do?

Shaft bed

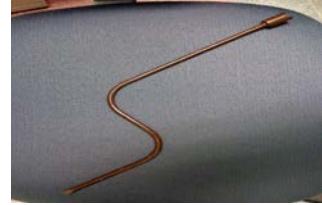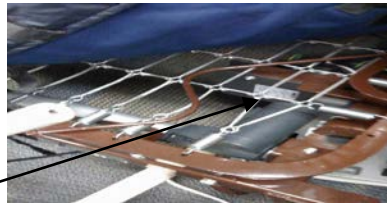

Locate Emergency Crank from under the mattress or under the foot end of the hospital bed.

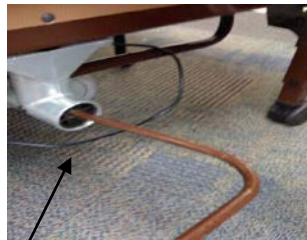


Place flat end of emergency crank in the center of the footboard to adjust bed height.



Use the other end of the crank for use with the motors.
To adjust the headspring and foot spring, select the desired motor end. The motor on the left side from the foot end will operate the head. The motor on the right will operate the foot.

Cord bed

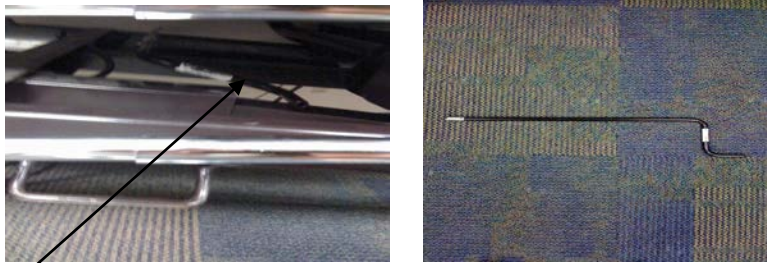


Locate the emergency crank from the center box of the bed.



Place the crank in the end of the motor. You will have 3 motors if it is a fully electric bed and 2 motors and a hand crank in the center if it is a semi-electric bed. Place the crank in the left motor to adjust the head and in the right motor to adjust the foot.

Bed Rails

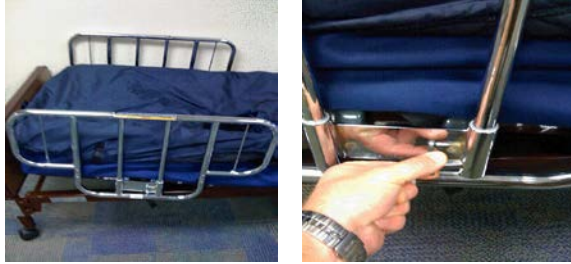

Half rails, split rails. Rails are only at the head end of the hospital bed.

Pull pin to raise or lower the side rail.



Full length hospital bed rails. Pull pin at head end and foot end to raise or lower rail.

Special Instructions