

## CADD Prizm Pump-TPN Profile

(Resetting the Reservoir Volume with a Bag Change)

### Set Up and Starting the Infusion

- Remove the battery from the pump, flip it over and put it back in the pump with the correct charges. (+ to + and - to -)




**\*\*\*\*\*Allow the pump to completely power up before attaching the IV tubing. You will hear 6 beeps to confirm the pump has successfully powered up. \*\*\*\*\***

### Spiking the Bag and Priming the Tubing

- Spike the IV tubing into the TPN bag port. Squeeze the bag to prime the tubing to remove the air until you see a drip on the end of the tubing. Allow the filter to hang straight down. You do not need to flip or rotate the filter.
- Remove the blue clamp from the cassette once the tubing is primed.

### Attaching the Cassette

- Attach the cassette hooks into the hinge pins on the pump. Place the pump on a hard surface and press the cassette down firmly.

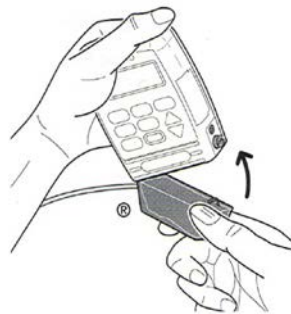

- Hold the pump securely. Using a coin, press in and turn the latch on the side of the pump counterclockwise until the mark on the latch is completely straight up.



- Screen will read **“High volume admin. set latched.”**

HIGH VOLUME ADmin  
set latched  
  
Next to continue

- Screen should read, **“STOPPED”**.
- Press **“NEXT”** to advance to the Reservoir Volume screen. The reservoir volume should be reset back to the amount of your TPN bag.
- If the “Res Vol” has not automatically reset, perform these steps:
  - Press **“ENTER”** while on the Reservoir Volume screen.
  - Screen will read, **“Reset Reservoir Volume to XXXX ml”?**
  - Press **Y** (Yes). Screen will read **“Reservoir Volume XXXX ml Resetting”. Infusion Profile has been reset.”**

Reset Reservoir  
Volume to \_\_\_\_\_ ml?  
  
Press Y or N

- **“Prime Tubing”?** Press **N** (No). **“Start the Pump?”** Press **N** (No).
- Press **“NEXT”** to advance through the program until you get to the **“Milliliters Given”** screen. Press **“ENTER”** to clear the amount infused/given back to **“0”** (zero).
- Press **Start**. Screen will read, **“Start the Pump?”** Press **Y** (Yes).

Start the Pump?  
  
Press Y or N

- Screen will read **“Starting Pump”** and then flip through each screen of the program.
- Screen will read **“TPN”** and the word **“RUNNING”** with an arrow displayed **↑** (for “ramp up”) or **→** (no “ramp up”). The “Res Vol” on the screen should read the amount of your bag/infusion.

### **Stopping the Infusion/Turning the Pump Off**

- Press **“STOP”**. Screen will read, **“Stop the Pump?”** Press **Y** or **N**.

Stop the Pump?  
  
Press Y or N

- Press **“Y** (Yes)”. Remove the bag and tubing from the pump and discard in the trash.
- Remove the battery from the pump, flip it over and put it back in the pump with the charges reversed. Lay the pump in a clean, safe place until you are ready to hook up the TPN.